1. Clean or prophylaxis with water and pumice powder;
2. Protection of soft tissues with vaseline or rubber dam;
3. Dry and apply the product during 2 - 3 minutes;
4. Wash.

SPECIAL CARES:

- Product is highly concentrated and caustic (pH 8,5). Protect soft tissues with vaseline or rubber dam.
- Presents secondary reaction regarding dental coloration (black stains) in the places where there are decays or structural failure, then it is not recommended to be used in permanent dentition.

- Protect yourself to avoid stains in the fingers and clothes when manipulating the product or its package.

NOTE: Repeat procedure at least once again and keep control over the patient.

SPECIAL RECOMMENDATIONS:

In case of accidental stains in skin apply on it a solution of Iodine tincture.

Dispose the product in accordance with local regulations, ensuring its complete distortion, preventing its reuse and environment damages.

CONTRAINDICATION:

CARIESTOP should not be used in patients is known to have sensibility to any of the ingredients listed.


WARNINGS:

Do not use material after expiry date.

The manufacturer is not liable for damages caused by improper use or use not provided at these instructions.

"KEEP MATERIAL OUT OF THE REACH OF CHILDREN".

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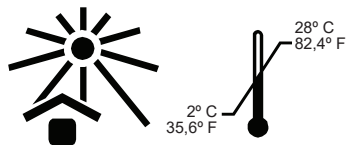ANVISA n o 10298550048.

COMPOSITION:

Fluoridric Acid, Silver Nitrate, Amonia Hydroxide and Deionized Water.

INDICATION:

Inhibiting agent from dental decay;

Prevention and inhibiting of dental decay in children;

For night feeding cases (when kids eat at night and do not brush their teeth appropriately);

Preventing secondary decays after amalgam restorations or melt down blocks;

Prevention of decays in pits and fissures;

As dentin desensitizer.

TECHNICAL INFORMATION:

Silver Fluoride solution 12%, 30% and 38% in amoniacal environment performance the strength of the enamel structure by forming CaF2 (calcium fluoride) and Ag3PO4 (silver phosphate), through its reaction with dental structure. Besides, there is the formation of proteical complexes of silver on the tooth surface. Preventing action is related to Fluoride ion and controlling action with Silver ion. Cariostatic action is obtained by the stimulation of dentin sclerosis, not allowing the progress of dental decay. Antimicrobial action is given mainly because its action over S. Mutans.

Researches and studies reported a reduction of 65% of decays in molars when cariostatics are applied every three months.

USE INSTRUCTIONS:

I - Intensive treatment:

4 applications: 1 per week during from 1 to 4 minutes according to the age.

II - Maintenance treatment:

1 application every 6 or 12 months.

Techniques or procedures